



**medi**

# medi FP10 / GP10 Quality compression garments on prescription

December 2010

medi. I feel better.





## medi

medi produce medical aids for a variety of conditions and life situations, but we don't see ourselves solely as a manufacturer. On the contrary, we handle our own research and development, technology, production, and marketing and sales in-house – always keeping our focus on the market, its demands and variables, and the medi vision.

## 2

medi

- 4 Contact
- 5 Quality Certificates

## 6

RAL quality hosiery on FP10/GP10

- 8 Overview
- 9 Indication for stockings
- 10 mediven
- 12 mediven plus
- 14 mediven elegance
- 16 mediven for men
- 18 mediven ulcer kit
- 20 mediven 95 armsleeves and gloves

## 22

British Standard (specification 40) Garments on FP10/GP10

- 22 duomed

## 24

Product overview charts

- 26 RAL quality
- 27 duomed



## We have the direct connection

Quality every time: in product, in service, in ordering and delivery.  
We are always here for you.  
This is your medi customer service team.

### Sales and customer service:

T 01432 373500  
F 01432 373510

You can visit our website: [www.mediuk.co.uk](http://www.mediuk.co.uk)  
or contact us by email: [fp10@mediuk.co.uk](mailto:fp10@mediuk.co.uk)  
or order by email: [sales@mediuk.co.uk](mailto:sales@mediuk.co.uk)

## Quality Certificates

medi's aim is to produce the best quality medical products in the world, guaranteeing effective solutions for a wide range of conditions.

### CE mark

Our medi medical products meet all of the requirements of the European Directive 93/42/EEC and of the German Medical Devices Act. In accordance with regulations, our products are subject to the conformity assessment procedures and are granted the CE mark. As a kind of “passport for industrial products”, the CE mark guarantees free trade within the European Economic Area (EEA).

### Quality label

Our medical compression garments are regularly checked by the Hohenstein Research Institute according to the quality and test provisions of the RAL-GZ 387. As confirmation that our products meet the quality conditions, you will find the quality label printed on the packaging.

### Textile trust quality seal

All compression stockings and garments are regularly examined and rated according to strict standards by independent test centres for any toxic substances potentially damaging to human health. Granting of the Oeko-Tex Standard 100 certificate guarantees: tested quality for the protection of your health.

### DIN EN ISO 9001:2000/DIN EN 13485

medi Bayreuth maintains the quality management system according to DIN EN ISO 9001: 2000 and DIN EN 13485. Since 1996, we have demonstrated our ISO 9000 status to our customers with a certificate. From this time, medi has exploited the opportunities offered by this standard to continually improve products from the business units Phlebology, Orthotics, Prosthetics and Hospital in order to organise safer and more efficient processes and procedures. We define customer satisfaction as a decisive measure of quality.

### Hohenstein testing device – HOSY

An accurate graduated compression profile is the most important feature of a medical compression garment.







## medi hosiery on prescription

Our comprehensive range of compression garments is now available on prescription.

We have supplied quality garments to Hospitals for many years – now we are ensuring that the best treatment can continue into the community by making an extensive product range available on FP10 / GP10.









## mediven®

The brand to remember for venous and oedema therapy.

medi thinks venous health right through to the end. Your patient can rely on mediven products; they are tailored for maximum compliance to guarantee the best possible treatment success.

We have the right compression garment for every indication. They are manufactured to the highest standard on modern machines using the best available yarns. All garments are clinically effective and are designed to be as comfortable as possible to encourage patient compliance. medi also supply a wide range of fitting aids and other accessories to assist with compression therapy.

medi creates first class brand quality with that certain something: the medi concept factor.



The tried and tested open toe RAL compression stocking.

Accurate graduated compression

To achieve the best possible clinical results mediven plus stockings are available in:

- 3 RAL compression classes
- 7 sizes
- 2 leg lengths
- Extra wide variants

Less made-to-measure

The wide range of sizes and variants available ensures that a vast majority of patients will fit into a standard garment without the need for made-to-measure options. Styles available include:

- Below knee
- Below knee extra wide
- Thigh length with silicone topband
- Thigh length extra wide with silicone topband
- Thigh with waist attachment

Mode of action

Tried and tested RAL quality compression, popular with both men and women. mediven plus provides reliable, accurate and comfortable compression.

Features











Patients will feel self confident wearing a stocking that not only provides accurate compression but is also soft, supple and hypo-allergenic.

The advanced construction of the yarn ensures mediven stockings are:

- Very comfortable
- Highly elastic
- Hard wearing and durable
- Suitable for washing and spin drying
- Latex free

medi stockings should be replaced after 6 months.

medi have a range of aids to assist in the application and removal of compression hosiery, please ask for details.

mediven plus								
Styles		AD	AD pet.	AD	AD pet.	AG/SN		
								
		calf	calf petite	calf x-wide	calf x-wide petite	thigh with topband		
Class 1	Code	DT 101	DT 102	DT 103	DT 104	DT 108		
Class 2	Code	DT 201	DT 202	DT 203	DT 204	DT 200		
Class 3	Code	DT 401	DT 402	DT 421	–	DT 400		
Styles		AG pet./SN	AG/SN	AG pet./SN	AG	AG pet.		
								
		thigh petite with topband	thigh x-wide with topband	thigh x-wide petite with topband	thigh with waist attachment	thigh petite with waist attachment		
				right		left	right	left
Class 1	Code	DT 106	DT 108X	DT 106X	DT 117	DT 119	DT 118	DT 120
Class 2	Code	DT 199	DT 200X	DT 199X	DT 207	DT 209	DT 208	DT 210
Class 3	Code	DT 399	DT 400X	DT 399X	DT 405	DT 406	DT 407	DT 408



Compression

RAL CCL 1, 2, 3

Product benefits

Pleasant on the skin, hard-wearing and opaque.



RAL class 1: 18 – 21 mmHg  
RAL class 2: 23 – 32 mmHg  
RAL class 3: 34 – 46 mmHg

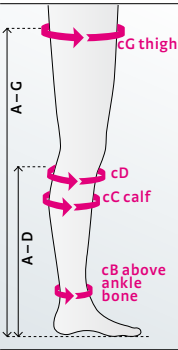
Extra wide thigh with topband available in beige only.

Colours  
mediven plus is available in beige and black on prescription.



length: size in cm		
	petite	regular
AD calf stocking	34 – 38	39 – 44
AG thigh stocking	62 – 71	72 – 83

Circumference: size in cms								
Size	I	II	III	IV	V	VI	VII	
cG thigh with silicone topband extra wide	49 – 57	53 – 62	57 – 67	61 – 72	65 – 77	69 – 81	73 – 85	
cG thigh with silicone topband	43 – 48	45 – 52	49 – 56	53 – 60	56 – 64	60 – 68	64 – 72	
cG thigh with waist attachment	43 – 57	45 – 62	49 – 67	53 – 72	56 – 77	60 – 81	64 – 85	
cC calf extra wide	34 – 38	37 – 41	40 – 44	42 – 46	44 – 48	46 – 50	48 – 52	
cC calf	28 – 34	30 – 37	33 – 40	35 – 43	37 – 46	39 – 49	41 – 51	
cB ankle	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32	



mediven®  
elegance

The soft and attractive closed toe option

The yarn and knit used to manufacture mediven elegance produces a sheer stocking that is clinically effective in RAL compression classes 1 and 2.

Available as:

- Below knee
- Thigh length with silicone topband
- Thigh length extra wide with silicone topband
- Tights

Mode of action

The sheer medical compression stocking with a unique soft and supple fit. The stocking offers an ideal combination of medical efficacy, wearing comfort and attractive appearance.

Features

Our compression stockings benefit from Clima-Comfort Technology. This brings your patients:

- Elasticity
- Softness
- Durability

mediven elegance sets new standards:

- The air-permeable leg section and the patent knit allow the skin to breath
- Hypo-allergenic materials minimise skin sensitivities
- Latex-free
- Softness: Uncoated inlaid threads and multi-filament fibres make them more comfortable to wear
- Closed toe with a highly elastic foot size

medi have a range of aids to assist in the application and removal of compression hosiery, please ask for details.



Compression

RAL CCL 1, 2

Product benefits

The elegant, fashionable stocking.











RAL class 1: 18 – 21 mmHg  
RAL class 2: 23 – 32 mmHg

Extra wide thigh with silicone topband available in beige only.


Colours  
mediven elegance is available in beige and black on prescription.



mediven elegance					
Styles		AD	AD pet.	AG/SN	AG pet./SN
					
		calf	calf petite	thigh with topband	thigh petite with topband
Class 1	Code	DT 185	DT 186	DT 189	DT 190
Class 2	Code	DT 285	DT 286	DT 289	DT 290
Styles		AG/SN	AG pet./SN	AT	AT pet.
					
		thigh extra wide with topband	thigh extra wide petite with topband	tights	tights petite
Class 1	Code	DT 189X	DT 190X	DT 191	DT 192
Class 2	Code	DT 289X	DT 290X	DT 291	DT 292

length: size in cm		
	petite	regular
AD calf stocking	34 – 38	39 – 44
AG full leg stocking	62 – 71	72 – 83
Closed toe stockings fit up to UK shoe size 9.		

Circumference: size in cms							
Size	I	II	III	IV	V	VI	VII
cG thigh with silicone topband extra wide	49 – 57	53 – 62	57 – 67	61 – 72	65 – 77	69 – 81	73 – 85
cG thigh with silicone topband	43 – 48	45 – 52	49 – 56	53 – 60	56 – 64	60 – 68	64 – 72
cG thigh for tights	43 – 57	45 – 62	49 – 67	53 – 72	56 – 77	60 – 81	64 – 85
cC calf	28 – 34	30 – 37	33 – 40	35 – 43	37 – 46	39 – 49	41 – 51
cB ankle	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32

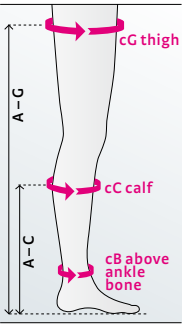



# mediven<sup>®</sup> for men

## Ribbed closed toe below knee stocking for the active patient

Venous disease affects men and women in roughly equal numbers - yet three times more women than men wear compression hosiery. To redress the balance medi have developed mediven for men, top quality compression garments that are designed to look, feel and be worn like conventional socks. And they still provide an accurate graduated compression profile. mediven for men – discreet and effective.

### Available in:

- RAL compression classes 1 and 2
- Calf length, standard and petite
- Closed toe for UK shoe sizes 3 - 11

### Sizes I to VII

### Mode of action

mediven for men is a below-knee stocking with a durable ribbed design that provides effective support whatever the activity. The construction of the high-tech fibres transports moisture to the outer surface of the stocking maintaining comfort even whilst active. The toe area is tapered to ensure a comfortable fit all day long. Latex-free.



### Compression

## RAL CCL 1, 2

### Product benefits



Effective support throughout all activities.



RAL class 1: 18 – 21 mmHg  
RAL class 2: 23 – 32 mmHg

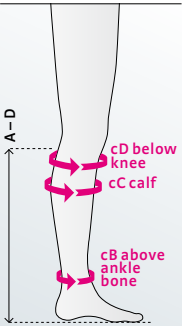
**Colours**  
mediven for men is available on prescription in the following colours.



mediven for men			
Styles		AD	AD pet.
			
		calf	calf petite
Class 1	Code	DT 342	DT 343
Class 2	Code	DT 340	DT 341
Length	2 lengths for AD		
Closed toe stockings fit up to UK shoe size 11.			

length: size in cms		
		petite      regular
AD calf stocking		33 – 38      39 – 44
Closed toe stockings fit up to UK shoe size 11		

Circumference: size in cms								
Size		I	II	III	IV	V	VI	VII
cD	below knee	27 – 33	29 – 36	32 – 39	34 – 42	36 – 45	38 – 48	40 – 50
cC	calf	28 – 34	30 – 37	33 – 40	35 – 43	37 – 46	39 – 49	41 – 51
cB	ankle	18 – 20	20 – 22	22 – 24	24 – 26	26 – 28	28 – 30	30 – 32





# mediven<sup>®</sup> ulcer kit

Cost effective treatment for crural venous ulcers.

The mediven ulcer kit helps you to manage your patients' venous ulcers without the bulk of a multi layer bandage system and the resources required to apply and maintain them. By using two stockings during the day, each providing 20mmHg, the clinically proven 40mmHg is applied at the ankle.

### How it works

The grey mediven ulcer liner stocking is applied over the top of any wound dressing. When active the patient wears the mediven plus stocking over the liner to give a total compression of 40mmHg at the ankle. The liner stockings contain interwoven silver which helps to reduce odour.

### mediven ulcer kit

- Ensures rapid application of accurate compression
- Allows patient to self manage pressure level
- Is less bulky than bandaging, improving comfort and mobility
- Offers an increase in patient confidence, independence and concordance
- Is machine washable

### Care instructions

mediven ulcer liner stocking can be washed at 95° C, mediven plus at 40° C. Both can be cool tumble dried.



### Product benefits

Unique product properties meet doctors' and patients' highest demands for reliable and safe graduated compression.



#### Inspection opening

The inspection opening under the foot of the mediven ulcer is primarily designed for patient safety. The circulation can be observed at any time and interdigital hygiene supported.



#### Soft end band

The soft end band gives a secure hold and does not constrict.



#### Instep/heel area

The knitted colour marking in the instep/heel area allows precise positioning in the heel and, together with the colour coded sizing band, allows rapid and reliable determination of stocking size. The shaped heel area provides optimum pressure reduction.

### Ordering information

The mediven ulcer kit can be ordered as a kit or as individual components.

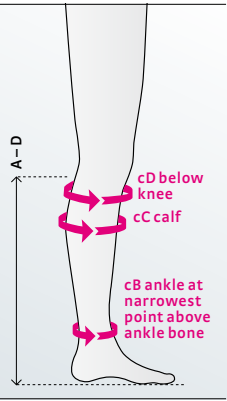
The kit contains:

- 2 x liners (understockings)
- 1 x mediven plus over stocking
- 1 x easy-on applicator

Additional mediven plus stockings can be ordered separately.

mediven ulcer kit				
Size	Leg measurements (circ. in cms)			Colour mark
	Ankle cB	Calf cC	Below knee cD	Heel
I	18 – 20	28 – 34	27 – 33	pink
II	20 – 22	30 – 37	29 – 36	blue
III	22 – 24	33 – 40	32 – 39	yellow
IV	24 – 26	35 – 43	34 – 42	red
V	26 – 28	37 – 46	36 – 45	grey
VI	28 – 30	39 – 49	38 – 48	green
VII	30 – 32	41 – 51	40 – 50	brown

length: Floor (A) to below knee (D) in cms		
	petite	regular
A-D calf stockings	34 – 38 cms	39 – 44 cms
ulcer kit	DT 276	DT 275
liner refill (2 per pack, universal length)	–	DT 277





# mediven® 95 armsleeve and gloves

## Quality upper limb compression

### Indications

For treatment in primary and secondary lymphoedema (ISL stages 1 & 2) of the arm with permanent or partial oedema of the hand. Post-operative and post-traumatic oedema and burns therapy.

### Mode of action

Growing knowledge encourages new developments. Since its market launch, mediven 95 has been continuously improved to benefit from new findings in lymphology.

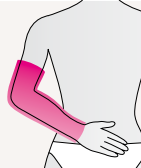

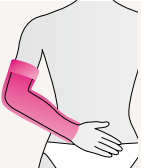
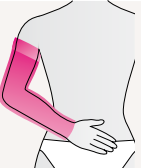





### Available in:

- Armsleeve
- Combined armsleeve
- Armsleeve with grip top
- Combined armsleeve with grip top
- Glove with compressive fingers

### Features

- Seamless, soft and attractive fibre
- Ease of handling, comfortable to wear
- Secure fit, non-slip knit
- Skin friendly, air permeable, fine-meshed
- Suitable for cool machine washing and cool tumble drying

Beige only  
Latex-free

mediven 95 armsleeve						
Styles		CG	CG long	CG/SN	CG/SN long	
						
		armsleeve	armsleeve long	armsleeve with grip top	armsleeve long with grip top	
CCL 1	Code	DT 714	DT 715	DT 746	DT 747	
CCL 2	Code	DT 716	DT 717	DT 750	DT 751	
Styles		AG	AG long	AG/SN	AG/SN long	AC <sub>1</sub>
						
		combined armsleeve	combined armsleeve long	combined armsleeve with grip top	combined armsleeve long with grip top	glove with compressive fingers
CCL 1	Code	DT 734	DT 735	DT 754	DT 755	DT 760
CCL 2	Code	DT 736	DT 737	DT 758	DT 759	DT 761



Compression

RAL CCL 1, 2

Product benefits

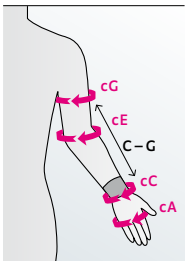


### Colours

mediven 95 is available in beige on prescription.



Standard sizes circumference in cms						
Sizes	II	III	IV	V	VI	VII
cG	22 – 29	25 – 32	29 – 36	33 – 40	37 – 44	41 – 48
cE	20 – 26	23 – 29	26 – 32	29 – 35	32 – 38	35 – 41
cC	13 – 15	15 – 17	17 – 19	19 – 21	21 – 23	23 – 25
cA	16 – 19	18 – 21	20 – 23	22 – 25	24 – 27	26 – 29
length C (wrist) to G (top of arm): standard 40 – 44 cm, long 44 – 48 cm						









duomed<sup>®</sup>  
British Standard  
(Specification 40)  
garments


The economic solution for venous problems.

Graduated compression elastic hosiery in compliance with British Standard (BS 6612:1985)

- Three compression classes
- Below knee and full thigh length available in all classes
- Machine washable
- Latex free
- Simple size selection
- Open toe in BS ccl 1, 2 and 3
- Closed toe in BS ccl 1 only

duomed open toe			
Styles		AD	AG
			
		calf	thigh
Class 1	Code	148	149
Class 2	Code	152	153
Class 3	Code	154	155
duomed closed toe			
Styles		AD	AG
			
		calf	thigh
Class 1	Code	150	151

duomed sizing chart				
Sizes	Mean ankle circumference	cB Ankle measurement	cC Calf measurement	cG Thigh measurement
1 (small)	20 cms	19 – 21 cms (7½ – 8½ ins)	28 – 34 cms (11 – 13½ ins)	42 – 56 cms (16½ – 22½ ins)
2 (medium)	23 cms	22 – 24 cms (8¾ – 9½ ins)	32 – 38 cms (12¾ – 15¼ ins)	48 – 62 cms (19¼ – 24¾ ins)
3 (large)	26 cms	25 – 27 cms (9¾ – 10¾ ins)	36 – 42 cms (14½ – 16¾ ins)	54 – 68 cms (21½ – 27¼ ins)
4 (x large)	29 cms	28 – 30 cms (11 – 12 ins)	40 – 46 cms (16 – 18½ ins)	60 – 74 cms (24 – 29½ ins)

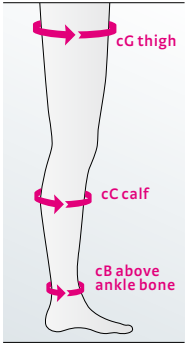

Compression

BS CCL 1, 2, 3

Product benefits

Economic wear for daily care.



Clima Fresh




CONFIDENCE  
IN TEXTILES

Tested for harmful substances  
according to Oeko-Tex Standard 100

BS class 1: 14 – 17 mmHg  
BS class 2: 18 – 24 mmHg  
BS class 3: 25 – 35 mmHg

Colours

Duomed is available in beige only.



Beige





## Product overview charts

medi understand the need for accurate graduated compression. To meet this demand we have a comprehensive range of garments available on prescription.

Stockings and armsleeves are designed to provide accurate graduated compression that reduces from the ankle (or wrist) towards the body. In order for the correct 'dose' of pressure to be applied it is important to select the correct compression class and then size of garment using the relevant measuring chart. An overview of the whole range follows but please remember that our customer services team are on hand during normal working hours if you need further assistance.

Simply telephone 01432 373500 during normal working hours and we promise our telephones will be answered by a person and not a machine!



product code overview

In order to prescribe the garment that your patient requires, state the mediven range (eg, plus, elegance etc), the compression class, style (eg, below knee), colour and size.

mediven plus – open toe, size I-VII							
Description		Below Knee	Below Knee Extra Wide	Thigh with Topband	Thigh Extra wide with Topband	Thigh with waist attachment	
						Right	Left
Class 1 18-21mmHg	Standard	DT101	DT103	DT108	DT108X*	DT117	DT119
	Petite	DT102	DT104	DT106	DT106X*	DT118	DT120
Class 2 23-32mmHg	Standard	DT201	DT203	DT200	DT200X*	DT207	DT209
	Petite	DT202	DT204	DT199	DT199X*	DT208	DT210
Class 3 34-46mmHg	Standard	DT401	DT421*	DT400	DT400X*	DT405	DT406
	Petite	DT402	N/A	DT399	DT399X*	DT407	DT408

\*Beige only

mediven elegance - closed toe, size I-VII					
Description		Below Knee	Thigh with Topband	Thigh Extra wide with Topband	Tights
Class 1 18-21 mmHg	Standard	DT185	DT189	DT189X*	DT191
	Petite	DT186	DT190	DT190X*	DT192
Class 2 23-32 mmHg	Standard	DT285	DT289	DT289X*	DT291
	Petite	DT286	DT290	DT290X*	DT292

\*Beige only

mediven for men - closed toe, size I-VII		
Description		Below Knee
Class 1 18-21 mmHg	Standard	DT342**
	Petite	DT343**
Class 2 23-32 mmHg	Standard	DT340**
	Petite	DT341**

\*\* Black or marine

mediven ulcer kit - size I-VII		
Description	Ulcer Kit (2x liner, 1x plus stocking)	Liner refill (2 per pack)
Standard	DT275	DT277
Petite	DT276	N/A

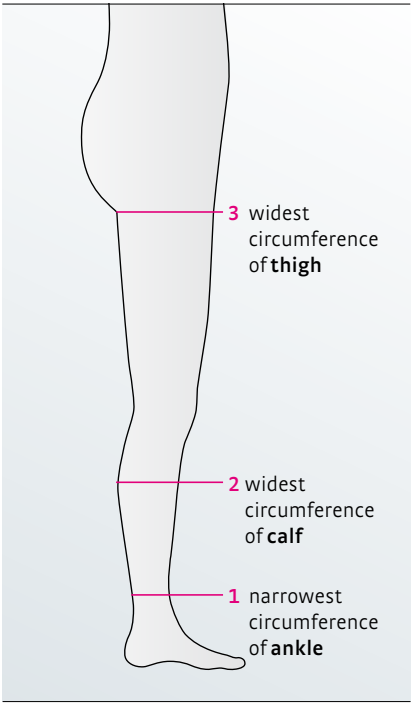
mediven 95 armsleeve and glove - sizes II-VII						
Description		Armsleeve	Armsleeve with griptop	Combined Armsleeve	Combined with griptop	Glove
Class 1	Standard	DT714	DT746	DT734	DT754	DT760
	Long	DT715	DT747	DT735	DT755	
Class 2	Standard	DT716	DT750	DT736	DT758	DT761
	Long	DT717	DT751	DT737	DT759	

Codes for ordering colours  
0 - Beige  
5 - Black  
3 - Marine (for men only)

Product style, code and sizing charts

Selection of correct size

- a) Measure the ankle, taking the circumference around the narrowest point (usually about 2 inches above the ankle bone).
- b) Select size of garment by referring to Column 1 in the chart below.
- c) For knee length stockings take one more measurement around the calf (the widest point below the knee) and check that this measurement falls within the range shown against the selected size in column 2 to ensure that the selected size will fit the patient.
- d) For thigh stockings take the calf measurement as referred to in (c) above and also a thigh measurement. Check that both these measurements fall within the ranges shown in columns 2 and 3 to ensure that the selected size will fit the patient.



duomed sizing chart				
		1	2	3
	Mean ankle circumference	Ankle circumference	Calf circumference	Thigh circumference
Small	20 cms	19 - 21 cms (7 - 8½ ins)	28 - 34 cms (11 - 13½ ins)	42 - 56 cms (16 - 22½ ins)
Medium	23 cms	22 - 24 cms (8¾ - 9½ ins)	32 - 38 cms (12¾ - 15¼ ins)	48 - 62 cms (19¼ - 24¾ ins)
Large	26 cms	25 - 27 cms (9¾ - 10¾ ins)	36 - 42 cms (14½ - 16¾ ins)	54 - 68 cms (21½ - 27¼ ins)
Extra Large	29 cms	28 - 30 cms (11 - 12 ins)	40 - 46 cms (16 - 18½ ins)	60 - 74 cms (24 - 29½ ins)

duomed style and code chart (Prescribable FP10)					
		Product codes / sizes			
duomed Class 1 (14 - 17 mmHg)	Light compression	Small	Medium	Large	Extra Large
Below knee	Open toe	148/1	148/2	148/3	148/4
Thigh	Open toe	149/1	149/2	149/3	149/4
Below knee	Closed toe	150/1	150/2	150/3	150/4
Thigh	Closed toe	151/1	151/2	151/3	151/4
duomed Class 2 (18 - 24 mmHg)	Medium compression	Small	Medium	Large	Extra Large
Below knee	Open toe	152/1	152/2	152/3	152/4
Thigh	Open toe	153/1	153/2	153/3	153/4
duomed Class 3 (25 - 35 mmHg)	Strong compression	Small	Medium	Large	Extra Large
Below knee	Open toe	154/1	154/2	154/3	154/4
Thigh	Open toe	155/1	155/2	155/3	155/4

medi UK Ltd.  
Plough Lane  
Hereford HR4 OEL  
Great Britain  
P 01432 373500  
F 01432 373510  
fp10@mediuk.co.uk  
www.mediuk.co.uk

